IDC Mikro






FRAMEWORK MANAGEMENT

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Deutsch	3 - 21
English	22 - 40

IDC

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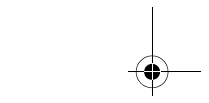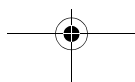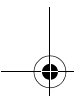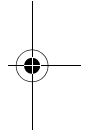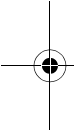DE

TABLE OF CONTENTS

- Translation of the original Instruction Manual -

EN

Table of Contents

1	Explanation of Symbols	23	8	Malfunctions, Repairs and Warranty	38
			8.1	Malfunctions	38
			8.2	Repairs	38
			8.3	Warranty	38
2	General Safety Instructions	24	9	Environmental Protection	38
3	Suitable Personnel	24	10	Technical Data and Accessories/Spare Parts	39
4	Machine Specifications	25			
4.1	Delivery Scope	25			
4.2	Intended Use	25			
4.3	CE Declaration of Conformity	25			
4.4	Components and Interfaces	26			
4.5	Accessories	27			
5	Installation	28			
5.1	Milling machine	28			
5.2	Maintenance unit	28			
5.3	Installation of the CAM software and milling machine	30			
5.4	Establishing the Network Connection	31			
5.5	Starting up the milling machine	31			
6	Application and Operation	33			
6.1	Changing the blank holder	33			
6.2	Function of the buttons	34			
6.3	Reference run	34			
6.4	Program execution	34			
6.5	Checking the workpiece during a running program	35			
6.6	Moving the the service position	35			
6.7	Moving to the change position for the blank holder	35			
7	Cleaning and Maintenance	36			
7.1	Cleaning	36			
7.2	Maintenance	36			

EXPLANATION OF SYMBOLS

1 Explanation of Symbols

Warning indications

Warning indications in the text are marked with a colour-backed triangle and boxed.

In case of hazards through electricity, the exclamation mark in the warning triangle is substituted by a lightning bolt.

Signal words at the beginning of a warning indication specify the type and severity of the consequences, if the measures to avert the hazard are not adhered to.

- **NOTE** means that property damage can occur.
- **CAUTION** means that light to fairly serious personal injury can occur.
- **WARNING** means that serious personal injury can occur.
- **DANGER** means that serious personal injury can occur.

Important information

Important information that do not lead to hazards for humans or property damage are marked with the icon aside and are also boxed.

Other symbols in the Manual

Symbol	Meaning
▷	Item of an operation description
–	Item of a list
•	Subitem of an operation description or a list
[3]	Numbers in square brackets refer to position numbers in graphics/figures

Tab. 1

Other symbols on the machine

Symbol	Meaning
	USB port (socket type B)
	Network connection (Ethernet)
	Slot for SD card
	Control cable connection, extraction
	Fuse
	Compressed-air connection

Tab. 2

EN

GENERAL SAFETY INSTRUCTIONS

EN

2 General Safety Instructions

When installing, starting-up and operating the machine, always observe the following safety instructions:

**CAUTION:**

A defective machine can lead to malfunctions!
On detection of damage or a functional defect of the machine:

- ▷ Label the machine as defective.
- ▷ Prevent further operation until the machine has been repaired.

**NOTE:**

Possible machine damage through insufficient extraction during dry processing!

- ▷ Operate the machine only with the extraction system recommended by the manufacturer or another compatible extraction system (see page 29).

**NOTE:**

- ▷ Switch the machine off when not in use or unsupervised for longer periods, e.g., overnight. This method of saving electrical energy also benefits the environment.

3 Suitable Personnel

**NOTE:**

Starting-up and operation of the machine may only be carried out by trained specialised personnel.

MACHINE SPECIFICATIONS

EN

4 Machine Specifications

4.1 Delivery Scope

- _ Milling machine
- _ Mains cable
- _ USB cable (5 m)
- _ Network cable (5 m)
- _ Control cable, extraction
- _ Compressed-air hoses (250 mm and 2000 mm)
- _ Cleaning hose with brush head
- _ Maintenance unit
- _ Service set for spindle
- _ CAM software package
- _ Test blank
- _ CNC cutter 2.5 (milling tool mounted in collet)
- _ Screw set
- _ Torque wrench, 4 Nm
- _ Blank holder 71
- _ Extractor adapter
- _ Slide caliper

The milling machine is delivered with a test blank inserted in the blank holder. Prior to delivery, a test object was factory-milled in order to check the zero point settings. Possibly remaining milled shavings/chips are not a reason for rejection. The test blank is provided for test milling and calibration of the milling machine.

- After unpacking, check the machine for completeness and possible transport damages. Please claim any transport damages immediately with your supplier.

4.2 Intended Use

The milling machine is a PC-controlled milling machine for the fabrication of dentures by means of:

- _ Dry processing of blanks made of
 - pre-sintered zirconium oxide
 - unsintered non-precious metals (NPMs)
 - wax
 - plastics

Using blanks and tools not approved by the manufacturer can damage the machine and make the product unusable. For such cases, the manufacturer shall assume no liability whatsoever.

Unauthorised modifications/alterations of the machine shall also void the warranty.

A compressed-air connection is required for operation of the milling machine (see page 29).

Dry processing is permitted only in conjunction with an appropriate extraction system (for technical requirements, see page 29).

4.3 CE Declaration of Conformity

In terms of design and performance, this product complies with the European Directives and the supplementary national requirements. Conformity has been confirmed with the CE marking.

EU Directives

2006/42/EG	2004/108/EG
------------	-------------

harmonized Standards

DIN EN ISO 12100	DIN EN 13128
------------------	--------------

DIN EN 61010-1	DIN EN 55011
----------------	--------------

DIN EN 61000-4-4

DIN EN 61000-4-5

Tab. 3

The product's Declaration of Conformity can be requested from the manufacturer.

MACHINE SPECIFICATIONS

EN

4.4 Components and Interfaces

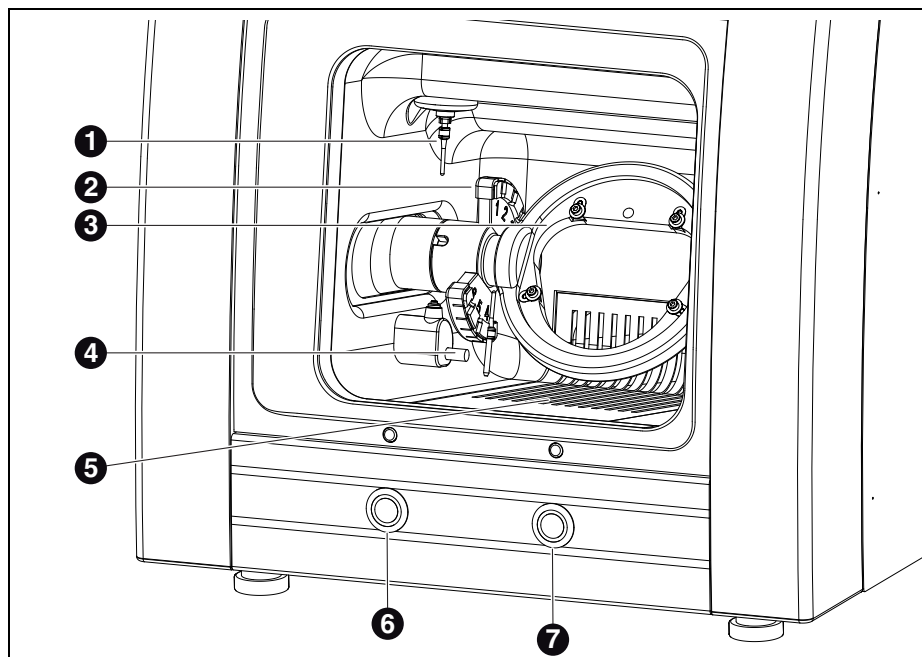
Fig. 1 Machine overview, front

- 1 Milling spindle
- 2 Tool holder
- 3 Blank holder
- 4 Measuring device, tool length
- 5 Collection tray/extraction duct
- 6 Operation button
- 7 Service button

MACHINE SPECIFICATIONS

EN

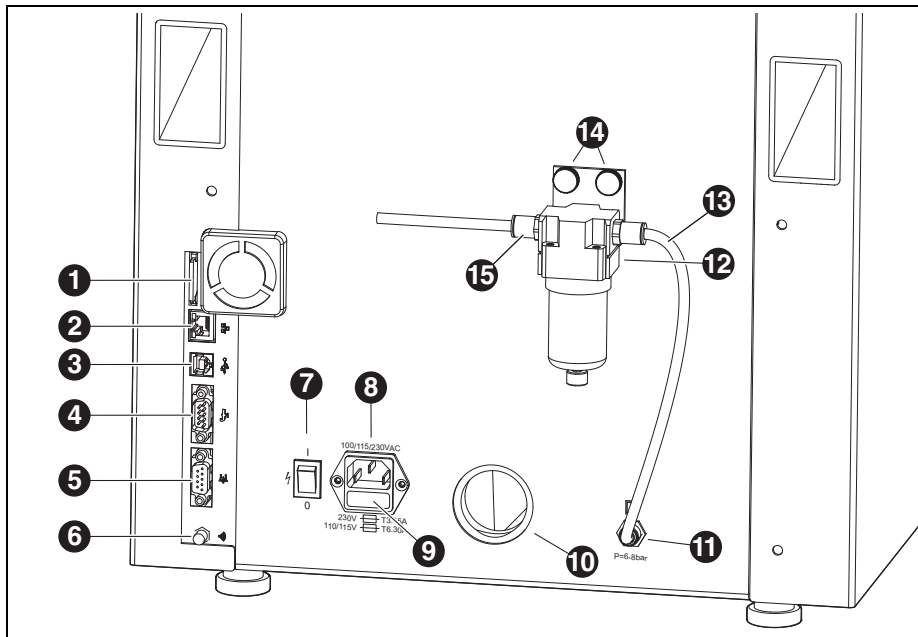
Fig. 2 Machine overview, rear

- 1 Slot for SD memory card
- 2 Ethernet interface
- 3 USB data port
- 4 Free connection
- 5 Control interface for extraction
- 6 WLAN antenna
- 7 Main switch
- 8 Power supply connection
- 9 Fuse compartment
- 10 Extraction-hose connection
- 11 Compressed-air connection, milling spindle
- 12 Maintenance unit
- 13 Compressed-air hose, 250 mm
- 14 Knurled screws
- 15 Compressed-air inlet

4.5 Accessories

Extraction

The extraction recommended by the manufacturer is required for dry processing with the milling machine. It is equipped with a special filter bag and a HEPA micro filter, which filters up to 99.97 % of fine dust particles (i.e. zirconium oxide dust particles) and is classified in the fine dust particle category M.

The separate socket for the control cable allows for automatic operation.

Milling tools

Various cutting tools are available for the milling machine (see catalogue).

INSTALLATION

EN

5 Installation

5.1 Milling machine

Countries with 100-115 V mains voltage

The machine is preset for a mains voltage of 230 V. When the mains voltage at the set-up location is 100-115 V:

- ▷ Pull out the fuse compartment ([9] in fig. 2 on page 27) below the mains plug.
- ▷ Replace both T3.15 A fuses with the provided T6.3 A fuses.

5.1.1 Set-up of the milling machine without supply unit

 Please refer to the separate milling machine set-up instructions.

- _ The machine is intended exclusively for use within dry, closed rooms.
- _ A clearance of at least 200 mm is to be observed sideways, to the rear and upwards. Minimum space requirements for the milling machine (incl. connections):
 - Width: 810 mm
 - Depth: 630 mm
 - Height: 740 mm
- _ The empty weight of the machine is 45 kg. The set-up surface must withstand an appropriate load.
- _ During operation, the room temperature must be between 18 °C and 30 °C; avoid large variations in temperature.
- ▷ Set up the machine on a sturdy work table or work bench (off the ground and clear of walls).
- ▷ Adjust the height of the machine by screwing the feet in or out, ensuring that the machine is positioned firmly on all four feet.

- ▷ Insert the extraction hose of the extraction system being used to the extraction-hose connection ([10] in fig. 2 on page 27) on the rear side of the milling machine.
- ▷ Connect the compressed-air hose on the rear side of the milling machine.

5.2 Maintenance unit

The spindle is equipped with a sealing-air feature. This sealing air prevents shavings/chips and dust from entering the spindle. The service unit filters possible contamination in the sealing air and thus prevents costly machine damage.

Installing the service unit

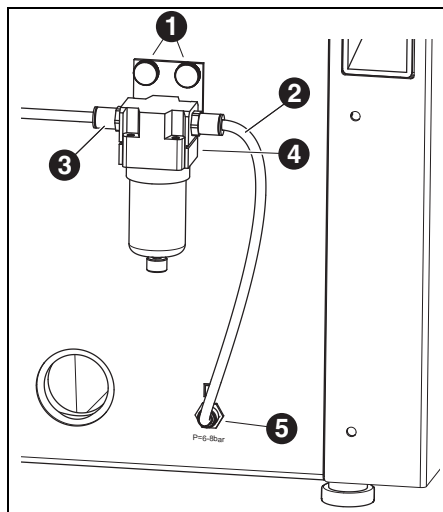
Fig. 3 Maintenance unit

- 1 Knurled screws
- 2 Compressed-air hose, 250 mm
- 3 Compressed-air inlet
- 4 Water separator
- 5 Compressed-air connection, milling spindle

INSTALLATION

EN

- ▷ Mount the service unit on the rear side of the machine using the knurled screws [1].
- ▷ Connect the short compressed-air hose [2] to the right-hand connection of the service unit and to the compressed-air connection [5] of the milling machine.
- ▷ Connect the long compressed-air hose to the compressed-air inlet [3] of the service unit and the other end to the compressed-air supply of the laboratory.

Requirements for the compressed-air supply

For operation, the compressed air must fulfil the following minimum requirements:

- _ Flow rate: 50 l/min (1.77 CFM)
- _ Pressure: 6 bar / 87 psi

Requirements on air purity:

- _ Solid contaminants: Class 3;
better than 5 µm for solids
- _ Water content: Class 4;
Max. pressure dew point +3°C
- _ Total oil content: Class 2;
Max. oil content 0.1 mg/m³

In continuous operation, the pressure and the air-flow rate are monitored by the milling machine. When the required values are fallen below, the milling program stops. Once the correct air values are available again, the milling program can be continued after fault acknowledgement in the software.

Requirements on air purity for the spindle

When another service unit than the one supplied is to be used, the following connection conditions according to DIN-ISO 8573-1 must be ensured:

- _ Solid contaminants: Class 3
(particle size max. 5 µm, particle density max. 5 mg/m³)
- _ Water content: Class 4
(max. pressure dew point +3 °C, water content max. 6000 mg/m³)
- _ Total oil content: Class 2;
(max. oil content 0.1 mg/m³)

Extraction requirements

When another extraction system is to be used, it must meet the following requirements:

- _ Suction capacity: 56.6 l/min
- _ Suitability for zirconium-oxide dusts
- _ HEPA microfilter (97.97 %), filter class H12, dust category M

The connection of the extraction hose to the machine is made with an adapter (in delivery scope). The adapter fits extraction hoses with a 38 mm interior diameter.

INSTALLATION

EN

5.3 Installation of the CAM software and milling machine

- ▷ Insert the data carrier into the drive of the PC.
The data carrier contains the installation file of the software.
- ▷ Open the file and follow the instructions in the installation.

Connecting the interfaces

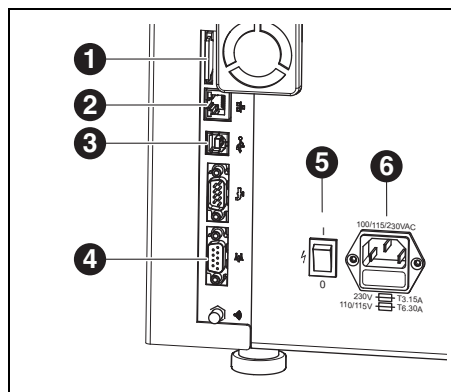
Fig. 4 Interfaces on the rear side of the milling machine

- 1 Slot for SD memory card
 - 2 Ethernet interface
 - 3 USB data port
 - 4 Control interface for extraction
 - 5 Main switch
 - 6 Fuse/power supply connection compartment
- ▷ Make sure that the main switch [5] of the milling machine is in the 0 position.
 - ▷ Check if the SD memory card is properly seated in the slot [1] on the rear side of the milling machine.
 - ▷ Connect the mains cable to power supply connection [6] of the milling machine as well as to a mains socket outlet.

- ▷ Connect the control cable to the control interface for extraction [4] on the rear side of the milling machine and to the interface on the rear side of the extraction system.

Connecting to the PC

- ▷ Connect a USB cable to the USB interface ([3] in fig. 4) of the milling machine and to a free USB port of the PC.
The maximum permitted length of the USB cable is 5 m.
- ▷ Insert the Dongle in a free USB port on the rear side of the PC.

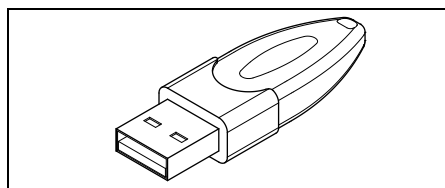
Fig. 5 CAM Dongle

INSTALLATION

5.4 Establishing the Network Connection

To establish the network connection, the milling machine and the PC must first be connected via USB.

- ▷ Connect the Ethernet interface ([2] in fig. 4) of the milling machine with the PC or with the network using a network cable.
The maximum permitted length of the network cable to the PC or to the network router is 100 m.
- ▷ Switch the milling machine on.
- ▷ Start the milling machine software.
- ▷ After initializing the milling machine in the software under *Settings > Communication*, set the connection type from USB to Ethernet and adapt the IP address, if required (see description of the milling machine software).
- ▷ Switch the milling machine off and then on again.
- ▷ Initialize the milling machine in the software.
- ▷ After the software has detected the Ethernet connection, remove the USB cable.

When operating via a network, up to eight milling machines can be connected and controlled via the CAM Software.
The PC and the milling machine must both be in the same sub-network.

5.5 Starting up the milling machine

The three standard cutters in tool stations 1-3 are mandatory for operation.

In addition to standard tools, sister tools can optionally be inserted in the tool holder. When a standard tool has reached its wear limit, the milling machine automatically accesses the respective sister tool.

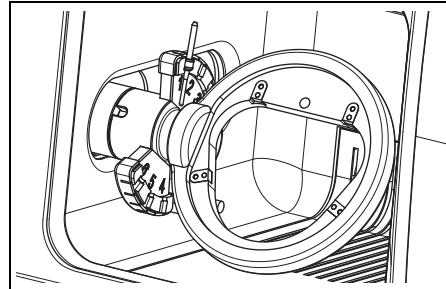
Fig. 6 Tool holder

- ▷ Insert the CNC cutter standard tools into the tool holder of the milling machine. Position the tool shank facing toward the outside.
 - CNC cutter 2.5 is in the spindle.
 - Set CNC cutter 1.0 to position 2.
 - Set CNC cutter 0.6 to position 3.
- ▷ Optionally: Insert the sister tools into the tool holder of the milling machine. Position the tool shank facing toward the outside.
 - Set CNC cutter 2.5 to position 4.
 - Set CNC cutter 1.0 to position 5.
 - Set CNC cutter 0.6 to position 6.

EN

INSTALLATION

EN

- ▷ Insert the selected blank into the holder in the milling machine. The high/thick side of the blank must face the inward.

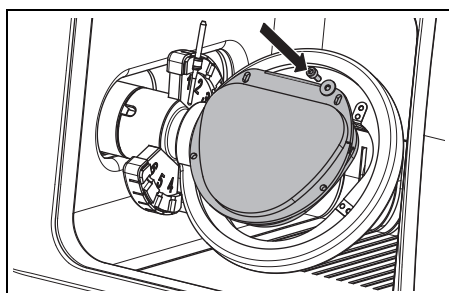
Fig. 7 Inserting a blank

- ▷ Fasten the blank using four screws and washers.

NOTE:

Do not damage the blank!

- ▷ Do not overtighten the screws! Avoid mechanical stress!

- ▷ Tighten the screws with the Allen key provided.

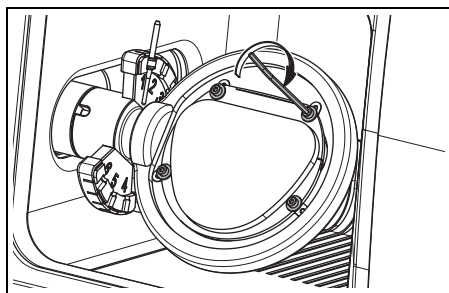
Fig. 8 Tightening a blank

- ▷ Shut the door of the milling machine.
- ▷ Switch the milling machine on via the main switch.
- ▷ Switch the PC on.

- ▷ Start the milling machine software.

When the software is started, a machine initialisation takes place.

If the milling machine is not detected:

- ▷ Restart the PC.

After re-positioning the milling machine, it is necessary to check the system axes.

- ▷ For this, observe the instructions for milling the test bodies.

APPLICATION AND OPERATION

EN

6 Application and Operation

WARNING:

Danger from shavings/chips or tool fragments being thrown from the machine!

- ▷ Make sure that the door of the milling machine is always closed during processing!

NOTE:

Insufficient milling results and increased wear when dry processing!

- ▷ Do not operate the milling machine without extraction!

- ▷ Pull off the blank holder toward the right.
- ▷ Clean the locating bore and surfaces.
- ▷ Check the pin of the blank holder and clean if required.
- ▷ Insert the blank holder with the pin to the stop in the locating bore.

NOTE:

Inaccurate working results or damage to the blank-holder fixture possible!

- ▷ Tighten fastening screw [1] with a torque of 4 Nm.

The torque wrench included in the delivery scope is appropriately set.

6.1 Changing the blank holder

In delivery condition, the milling machine is equipped with a 71 blank holder.

To change the blank holder:

- ▷ With the door closed, press operation button [A] for at least three seconds.

The current tool is placed down. The blank holder [2] moves to the change position.

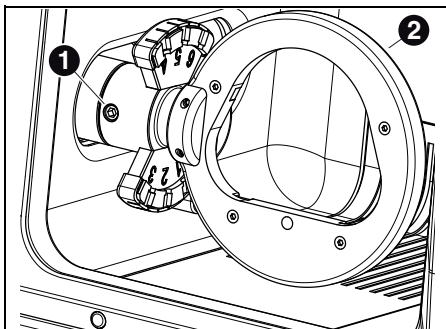
Fig. 9

1 Fastening screw for blank holder

2 Blank holder

- ▷ Loosen fastening screw [1].

- ▷ Tighten fastening screw [1] with the supplied torque wrench.
- ▷ Briefly press operation button [A].
The machine advances to the starting position.

APPLICATION AND OPERATION

EN

6.2 Function of the buttons

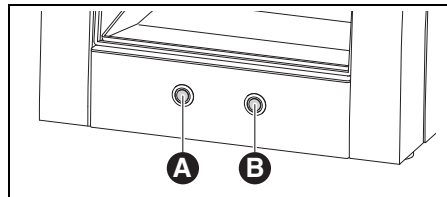
Fig. 10 Buttons

Key	Door shut	Door open
Operation button [A]	Starts or stops the reference run; starts or stops the current program; activates or deactivates the control position; activates the change position for the blank holder	Activates or deactivates the extraction system
Service button [B]	Interrupts the current program; activates the service position; activates the tool-changing position; starts a cleaning run of the inspection window	Activates or deactivates the sealing air

Tab. 4

Instead of with the buttons, these functions can also be actuated via the software.

6.3 Reference run

After switching on, a reference run must be carried out to determine the mechanical zero points. Operation button [A] flashes.

To start the reference run:

- ▷ Briefly push operation button [A].

The reference run is started. Operation button [A] lights up.

The reference run can be stopped:

- ▷ Briefly push operation button [A].

The reference run is stopped. Operation button [A] flashes.

To continue the reference run:

- ▷ Briefly push operation button [A].

The reference run is continued at the position where it was stopped. Operation button [A] lights up.

On completion of the reference run, the machine returns to the starting position. Operation button [A] goes out.

6.4 Program execution

Before executing a program, it must first be loaded (see CAM Software).

- ▷ Briefly push operation button [A].

The program is started. Operation button [A] lights up.

The program can be stopped:

- ▷ Briefly press operation button [A] or service button [B].

-or-

- ▷ Open the door.

The program is stopped. Operation button [A] flashes.

APPLICATION AND OPERATION

EN

To continue the program:

- ▷ With the door closed, briefly press operation button [A].
The program is continued at the position where it was stopped. Operation button [A] lights up.

On completion of the program, the machine returns to the starting position. Operation button [A] goes out.

6.5 Checking the workpiece during a running program

- ▷ Push and hold operation button [A] for at least three seconds.
The program is interrupted. The spindle returns to the rear, the current tool is placed down and the workpiece is turned around. While travelling, operation button [A] is lit; after reaching the position, it flashes.
- ▷ Visually check the workpiece.
- ▷ Briefly push operation button [A].
The tool is picked up again. The spindle travels to the position, where the program was interrupted. While travelling, operation button [A] is lit; after reaching the position, it flashes.
- ▷ Briefly push operation button [A].
The program is continued. Operation button [A] lights up.

6.6 Moving the the service position

The service position can be moved to from the starting position or from an interrupted program.

- ▷ Briefly press service button [B].
The current tool is placed down and the spindle returns to the service position. The spindle is released; the motor is blocked. While travelling, operation button [A] is lit; after reaching the position, it flashes.

To exit the service position:

- ▷ Briefly press operation button [A].
The tool is picked up again. The spindle travels back to the starting position or to the position, where the program was interrupted. While travelling, operation button [A] is lit; after reaching the position, it flashes.

6.7 Moving to the change position for the blank holder

The change position for the blank holder can only be moved to from the starting position.

- ▷ Push and hold operation button [A] for at least three seconds.
The current tool is placed down. The blank holder is turned in such a manner that the fastening screw becomes accessible. While travelling, operation button [A] is lit; after reaching the position, it flashes.

To exit the change position for the blank holder:

- ▷ Briefly press operation button [A].
The machine moves to the starting position.

CLEANING AND MAINTENANCE

EN

7 Cleaning and Maintenance

7.1 Cleaning

The milling machine must be cleaned after each working day.

NOTE:

Zirconia dust is abrasive!

- ▷ Clean the window in the door carefully, using only a soft cloth; do not scrub or scour.

- ▷ Open the door of the milling machine.
- ▷ Press operation button [A] on the milling machine.
The extraction system is activated. The operation button lights up.
- ▷ Press service button [B] on the milling machine.
The sealing air at the spindle is activated. The service button lights up.
- ▷ Remove the collecting plate.
- ▷ Vacuum out the milling chamber.
- ▷ Briefly press service button [B].
The sealing air is deactivated.
- ▷ Briefly press operation button [A].
The extraction system is deactivated.

7.2 Maintenance

7.2.1 Weekly maintenance

NOTE:

Danger of damaging the machine!

- ▷ Do not use compressed air, ultrasonic or a steam jet for the cleaning!

Prior to maintenance:

Before maintaining the spindle, the spindle sealing air must be activated:

- ▷ Briefly press service button [B].
The current tool is placed down and the spindle returns to the service position. The spindle is released; the motor is blocked.
- ▷ Open the door of the milling machine.
- ▷ Press service button [B] on the milling machine.
The sealing air is activated. The service button lights up.

After the maintenance

When the maintenance is completed:

- ▷ Shut the door of the milling machine.
- ▷ Briefly push operation button [A].
The tool is picked up again. The spindle travels back to the starting position or to the position, where the program was interrupted.

Spindle

The provided spindle service set and a milling tool of your choice are required for maintenance of the spindle.

In order to clean the spindle:

- ▷ Remove any dust from the cap and the spindle using a dry, soft cloth. Do not use cleaning agents!

CLEANING AND MAINTENANCE

To remove the collet:

- ▷ Apply the spindle wrench to the collet.

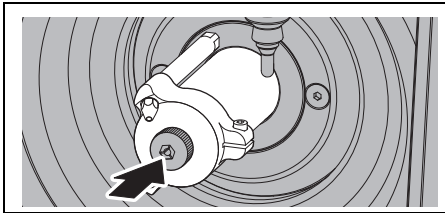
Fig. 11

- ▷ Insert a milling tool into the open collet.
- ▷ Using the spindle wrench, unscrew the collet out of the spindle.

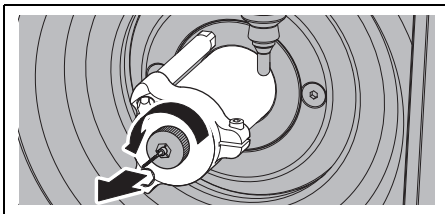
Fig. 12

- ▷ Clean the collet holder.
The collet holder must be free of shavings/chips and debris.

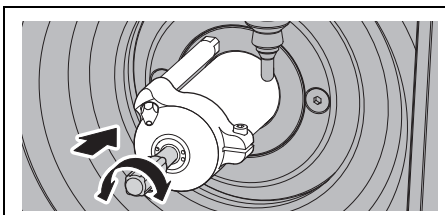
Fig. 13

- ▷ Clean the inside of the collet with the fine collet brush.

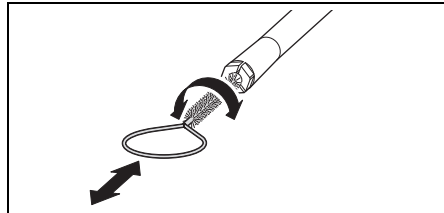
Fig. 14

NOTE:

Only use the provided collet grease when maintaining the collet cone.

- ▷ Apply a slight coat of collet grease to the threads of the collet.
- ▷ Screw the cleaned and lubricated collet in reverse order to the stop back into the spindle.

CAUTION:

Danger of tool breakage!

- ▷ Make sure to remove the milling tool from the collet chuck!

When the milling tool was from the milling machine tool holder:

- ▷ Place the tool properly positioned back into the tool holder.

EN

MALFUNCTIONS, REPAIRS AND WARRANTY

EN

Maintenance unit

- ▷ Visually check the service unit for:
 - Visible particles, debris or particulate matter within the acrylic glass (signs for contaminated compressed air)
 - Amber discolouration of the filter element (sign of oil in the compressed air)
 - Water in the filter (sign of water in the compressed air)

Once any of the listed impairments is given:

- ▷ Correct the problem with the compressed air.
- ▷ Replace the complete service unit.

7.2.2 Monthly Maintenance

Checking the zero point

It is recommended to check the zero point once every month.

- ▷ Mill a test cube from the test blank.

The test blank is to be used exclusively for test milling and calibration of the milling machine.

7.2.3 External maintenance

External maintenance is necessary after a defined number operating hours. The milling machine software will inform you appropriately.

8 Malfunctions, Repairs and Warranty

8.1 Malfunctions

In case of malfunctions:

- ▷ Restart the software.
- ▷ Restart the milling machine.
- ▷ Restart the PC.

8.2 Repairs

Repairs may be carried out only by trained, specialised personnel.

8.3 Warranty

The warranty complies with the statutory provisions. For more information, please refer to our General Terms of Business.

9 Environmental Protection

Packaging

In terms of packaging, the manufacturer participates in country-specific recycling systems, which ensure optimal recycling.

All packaging materials used are environmentally friendly and recyclable.

Waste equipment

Waste equipment contain materials that should be sorted for environmental-friendly recycling.

- ▷ After the service life has elapsed, dispose of the machine via the public disposal systems in an environmentally-friendly manner.

The plastic parts are labelled. Thus, they can be sorted for recycling or disposal.

TECHNICAL DATA AND ACCESSORIES/SPARE PARTS

10 Technical Data and Accessories/Spare Parts


Subject to changes.

EN

Technical data - milling machine

	Unit	Value
Art. No.	–	179300-IDC
Dimensions (D × W × H)	mm	588 × 516 × 755
Weight	kg	45
Electrical connections	V/Hz	100-230/50-60
Power output	W	250
Fuse (time-delayed)	A	3.15 / 6.3
Axes	–	4
Accuracy	μm	< 10
Allowable temperature range (operation)	°C	18 - 30
Spindle torque	Ncm	4
Spindle speed	rpm	60,000
Collet diameter	mm	3
Sound pressure level	dB(A)	60
Interfaces	–	USB / Ethernet / slot for SD card / Extraction control system
Compressed-air connection	–	Dry, clean compressed air
_ Pressure	bar	6
_ Flow rate	l/min	50
Air purity:		
_ Solid contaminants	–	Class 3; better than 5 μm for solids
_ Water content	–	Class 4; Max. pressure dew point +3 °C
_ Total oil content:	–	Class 2; max. oil content 0.1 mg/m ³
Extraction	–	Suitable for zirconium oxide dust particles, HEPA microfilter (97.97 %), filter class H12, dust category M
_ Diameter of extraction connection	mm	38
_ Suction capacity	l/min	56.6

Tab. 5

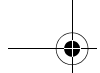TECHNICAL DATA AND ACCESSORIES/SPARE PARTS













EN





Accessories/Spare parts






Art. No.	Designation
179210	Maintenance unit
179211	Screw set for blank holder
179216	Threaded inserts for blank holder
179251	Cover, tool holder
179257	Extraction hose with brush
1156901	Fuse (T 3.15 A / 250 V)
179186	Blank holder 71

Tab. 6

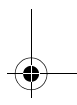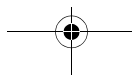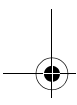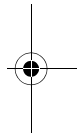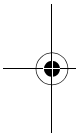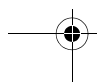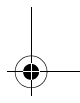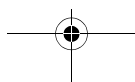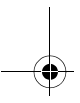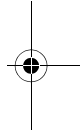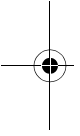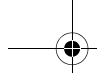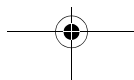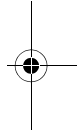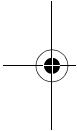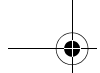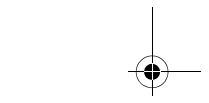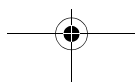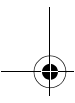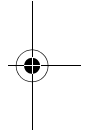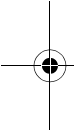32919-FB-IDC 2015-04-27

Made in the European Union

CE

Manufacturer | Hersteller
Distribution | Vertrieb

Amann Girrbach AG
Herrschaftswiesen 1
6842 Koblach | Austria
Fon +43 5523 62333-105
Fax +43 5523 62333-5119
www.idc.dental

IDC