Power Tower View



Installation and operating instructions





Contents

					6.2 6.3	Hose materials	15
lm	norta	nt information				and pipes	15
	-	nt information	0		6.4	Information about electrical con-	
1		this document	3			nections	15
	1.1	Warnings and symbols	3		6.5	Information about connecting	40
	1.2	Copyright information	4			cables	16
2		/	4	7		llation	16
	2.1	Intended purpose	4		7.1	Combining devices safely	16
	2.2	Intended use	4		7.2	Setting up the suction module	16
	2.3	Improper use	4		7.3	Install the compressed air mod-	
	2.4	General safety information	4			ule	17
	2.5	Specialist personnel	5		7.4	Installing the compressor unit	17
	2.6	Electrical safety	5		7.5	Positioning the unit	17
	2.7	Only use original parts	5		7.6	Installing a second compressor	40
	2.8	Transport	5			unit (optional)	19
	2.9	Disposal	5		7.7	Install the drying unit	19
					7.8	Additional work for 230 V variant .	21
Œ	<u> </u>				7.9	Install the amalgam separator	00
					7.40	(optional)	23
Pr	oduct	t description			7.10	Connecting the unit	23
3	Overv	iew	6		7.11	Two devices in a single com-	24
	3.1	Scope of delivery	6			pressed air network	
	3.2	Optional items	6	8		missioning	25
	3.3	Wear parts and replacement			8.1	Monitoring the unit via the net-	
		parts	6			work	25
4	Techn	ical data	8				
	4.1	400 V 3~ variant	9		5		
	4.2	230 V 1~ variant	9	He	age		
	4.3	Type plate	11		•		
	4.4	Declaration of conformity	12	9	•	ating the touch screen	27
_		•			9.1	Navigating	27
5	•	N/at a satisfa a satisfa	13		9.2	Using menus	27
	5.1	Wet suction system	13		9.3	Calling up messages on the	07
	5.2	Dry suction system	13			touch screen	27
	5.3	Functional description	14	10	Opera	ation	28
	5.4	Touch display	14		10.1	Switch on the unit	28
					10.2	Device information	28
					10.3	Configuring the units	28
۸۰		N.			10.4	Amalgam separator	28
	semb	•		11	Maint	tenance	29
6	-	rements	15	•••	11.1	Compressor unit and drying unit.	29
	6.1	Pipe materials	15			compression and anymig unit.	20

?

Troubleshooting

12	Tips for operators and service techni-	
	cians	31



Appendix

Αþ	pena	IIX	
13	Conne	ection media plan	32
14	Electr	ical connection plan	33
	14.1	Version 400 V, 3N~	33
	14.2	Version 230 V, 1N~	34
15	Menu	structure	35
	15.1	Settings	35

Important information

1 About this document

These installation and operating instructions represent part of the unit.



If the instructions and information in these installation and operating instructions are not followed, Dürr Dental will not be able to offer any warranty or assume any liability for the safe operation and the safe functioning of the unit.

The German version of the installation and operating instructions is the original manual. All other languages are translation of the original manual. These operating instructions are valid for Power Tower View, order number: A949444000, A949443000, A949453000, A949242000, A949252000 and A949243000.

1.1 Warnings and symbols

Warnings

The warnings in this document are intended to draw your attention to possible injury to persons or damage to machinery.

The following warning symbols are used:



General warning symbol



Warning - dangerous high voltage

The warnings are structured as follows:



SIGNAL WORD

Description of the type and source of

Here you will find the possible consequences of ignoring the warning

Follow these measures to avoid the danger.

The signal word differentiates between four levels of danger:

- DANGER

Immediate danger of severe injury or death

WARNING

Possible danger of severe injury or death

- CALITION

Risk of minor injuries

- NOTICE

Risk of extensive material/property damage

Other symbols

These symbols are used in the document and on or in the unit:



Note, e.g. specific instructions regarding efficient and cost-effective use of the unit.



Observe the operating instructions.



CE labelling



Manufacturer



Serial number



Order number



Dispose of correctly in accordance with EU Directive 2012/19/EU (WEEE).



Dispose of the packaging material in an environmentally responsible manner.



Wear protective gloves.



Wear protective goggles.



3N Three phase alternating current with neutral conductor

1N Single phase alternating current with neutral conductor



Disconnect all power from the unit.



Monitor ambient conditions



1.2 Copyright information

All circuits, processes, names, software programs and units mentioned in this document are protected by copyright.

The Installation and Operating Instructions must not be copied or reprinted, neither in full nor in part, without written authorisation from Dürr Dental

2 Safety

Dürr Dental has designed and constructed this unit so that when used properly and for the intended purpose it does not pose any danger to people or property.

Despite this, the following residual risks can remain:

- Personal injury due to incorrect use/misuse
- Personal injury due to mechanical effects
- Personal injury due to electric shock
- Personal injury due to radiation
- Personal injury due to fire
- Personal injury due to thermal effects on skin
- Personal injury due to lack of hygiene, e.g. infection

2.1 Intended purpose

Noise-insulated central supply of dental compressed air as well as negative pressure for dental suction. Optionally, water and amalgam can be separated.

2.2 Intended use

Central installation of the practice supply units in any room of a dental practice. Due to the noise insulation, the room can also be used for other purposes.

2.3 Improper use

Any use of this appliance / these appliances above and beyond that described in the Installation and Operating Instructions is deemed to be incorrect usage. The manufacturer cannot be held liable for any damage resulting from incorrect usage. The operator will be held liable and bears all risks

2.4 General safety information

- Always comply with the specifications of all guidelines, laws, and other rules and regulations applicable at the site of operation for the operation of this unit.
- Check the function and condition of the unit prior to every use.
- > Do not convert or modify the unit.
- Comply with the specifications of the Installation and Operating Instructions.
- The Installation and Operating Instructions must be accessible to all operators of the unit at all times.

2.5 Specialist personnel

Operation

Unit operating personnel must ensure safe and correct handling based on their training and knowledge.

Instruct or have every user instructed in handling the unit.

The following groups are not permitted to operate or use a commercially operated unit:

- People without the necessary experience and knowledge
- People with reduced physical, sensory or mental capabilities
- Children

Installation and repairs

Installation, readjustments, alterations, upgrades and repairs must be carried out by Dürr Dental or by qualified personnel specifically approved and authorized by Dürr Dental.

2.6 Electrical safety

- Comply with all the relevant electrical safety regulations when working on the unit.
- Never touch the patient and unshielded plug connections on the device at the same time.
- > Replace any damaged cables or plugs immediately.

2.7 Only use original parts

- Only use accessories and optional items that have been recommended or specifically approved by Dürr Dental.
- Only use only original wear parts and replacement parts.

2.8 Transport

The original packaging provides optimum protection for the unit during transport.

If required, original packaging for the unit can be ordered from Dürr Dental.



Dürr Dental will not accept any responsibility or liability for damage occurring during transport due to the use of incorrect packaging, even where the unit is still under guarantee.

- > Only transport the unit in its original packaging.
- > Keep the packing materials out of the reach of children.

2.9 Disposal



The unit must be disposed of properly. Within the European Union, the unit must be disposed of in accordance with EU Directive 2012/19/EU (WEEE).

If you have any questions about the correct disposal of parts, please contact your dental trade supplier.



The unit may be contaminated. Instruct the company disposing of the waste to take the relevant safety precautions.

Decontaminate potentially contaminated parts before disposing of them.



An overview of the waste keys for Dürr Dental products can be found in the download area at:

www.duerrdental.com
Document no.: P007100155



Product description

3 Overview

3.1 Scope of delivery

The following items are included in the scope of delivery (possible variations due to country-specific requirements and/or import regulations):

- Basic unit
- 400 V electronics box
- 400 V VS 1200 S
- 400 V Duo unit
- 20 Loressure tank
- Universal membrane dryer
- OroCup Care system

- Basic unit
- 400 V electronics box
- 400 V VS 900 S
- 400 V Duo unit
- 20 I pressure tank
- Universal membrane drver
- OroCup Care system

- Basic unit
- 400 V electronics box
- V 1200 S 400 V
- 400 V Duo unit
- 20 I pressure tank
- Universal membrane dryer

- Basic unit
- 400 V electronics box
- V 900 S 400 V
- 400 V Duo unit
- 20 I pressure tank
- Universal membrane dryer

- Basic unit
- 230 V electronics box
- 230 V/VS 600
- 230 V Duo unit
- 20 I pressure tank
- Universal membrane drver
- OroCup Care system

- Basic unit
- 230 V electronics box
- V 600 230 V
- 230 V Duo unit
- 20 l pressure tank
- Universal membrane dryer

- Basic unit
- 230 V electronics box
- 230 V VS 900
- 230 V Duo unit
- 20 I pressure tank
- Universal membrane dryer
- OroCup Care system

3.2 Optional items

The following optional items can be used with the device:

 Pressure reducer
 6040-992-00

 Amalgam separator
 7805-100-50

 Compressor unit 400 V 3~
 A949300059

 Compressor unit 230 V 1~
 A949300060

 Retrofit unit 230 V including valve block
 A949100098

 Valve block for 230 V 1~
 A949300092

3.3 Wear parts and replacement parts

The following working parts must be replaced at regular intervals (refer also to "Maintenance"); these articles do not bear the CE mark:

 Air intake filter
 0832-982-00

 Fine filter
 1610-121-00

 Virus bacteria filter
 1650100172

 Sintered filter
 1650-101-00

 Coalescence filter
 1650200323



To configure the required filters or filter sets, you can also use our filter configurator at:

www.duerrdental.com/filterkonfigurator



Any repairs exceeding routine maintenance may only be carried out by qualified personnel or our service.



Information about replacement parts is available from the portal for authorised specialist dealers at:

www.duerrdental.net



Hz	50
	I
	IP 20
	S1
mm	640 x 2080 x 610
I	20
°C	<5
MPa	1
MPa	0.5
MPa	0.9
erature of +40°C	
mm	7.2
mm	Ø 50
mm	Ø 50
mm	Ø 20
	Ethernet
	IEEE 802.3u
Mbit/s	100
	RJ45
	Auto MDI-X
	≥ CAT5
	RJ45
	≥ CAT5
ansport	
°C	-25 to +55
	mm I °C MPa MPa MPa MPa erature of +40°C mm mm mm mm mm Mbit/s



. 10 to . 10
+10 to +40
70
< 2000

4.1 400 V 3~ variant

Delivery at 5 bar (0.5 MPa)*	l/min	123
Noise level **	dB(A)	55
Noise level ***	dB(A)	53
Current reduction without CA 4 amalgam separator	Α	-0.2
Current reduction without second Duo unit	А	-2.5

^{* 1} compressor unit without drying unit

PTV A949444000 and A949443000 with second Duo unit and CA 4

General technical data		
Weight	kg	275
Electrical data		
Rated voltage	V	400, 3N~
Nominal current at 10 bar (1.0 MPa)	А	9.0

PTV A949454000 and A949453000 with second Duo unit

General technical data			
Weight	kg	265	

Electrical data		
Rated voltage	V	400, 3N~
Nominal current at 10 bar (1.0 MPa)	А	8.8

4.2 230 V 1~ variant

Delivery at 5 bar (0.5 MPa)*	l/min	115
Noise level **	dB(A)	55
Noise level ***	dB(A)	55
Current reduction without CA 4 amalgam separator	А	-0.5
Current reduction without second Duo unit	А	-6
Max. permitted mains impedance in accordance with IEC 61000-3-11	Ω	0.24

^{**} At an ambient temperature of +40 °C, noise level in accordance with ISO 3744

At an ambient temperature of +35 °C, noise level in accordance with ISO 3744



- 1 compressor unit without drying unit
- * At an ambient temperature of +40 °C, noise level in accordance with ISO 3744
- *** At an ambient temperature of +35 °C, noise level in accordance with ISO 3744

PTV A949242000 with second Duo unit

75
, 1N~
8
1

PTV A949252000 with second Duo unit


General technical data		
Weight	kg	265
Electrical data		
Rated voltage	V	230, 1N~
Nominal current at 9 bar (0.9 MPa)	А	17.5

PTV A949243000 with second Duo unit

General technical data			
Weight	kg	275	
Electrical data			
Rated voltage	V	230, 1N~	
Nominal current at 9 bar (0.9 MPa)	Α	20.4	

4.3 Type plate

The type plate is located at the bottom right of the housing.

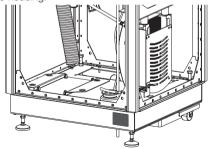

ΕN

4.4 Declaration of conformity

Name of manufacturer: DÜRR DENTAL SE
Address of manufacturer: Höpfigheimer Str. 17

74321 Bietigheim-Bissingen

Germany

Name of product: Power Tower View

We hereby declare that the above product meets all applicable requirements of the directives listed below:

- Machinery directive 2006/42/EC in its current version.

Name of the person with responsibility for compiling all technical documentation:

J. Gramann (Head of Research and Development).

The address is identical to the address of the manufacturer.


- Low voltage guideline 2014/35/EU in its valid version.
- Electromagnetic Compatibility (EMC) Directive 2014/30/EU in its current version.
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment in its current version.

Power Tower View	REF A949444000 REF A949443000 REF A949454000 REF A9499453000 REF A949242000
	REF A949252000 REF A949243000


J. Gramann on behalf of Dürr Dental Head of Research and Development O. Lange on behalf of Duerr Dental Head of Quality Management

5 Operation

5.1 Wet suction system

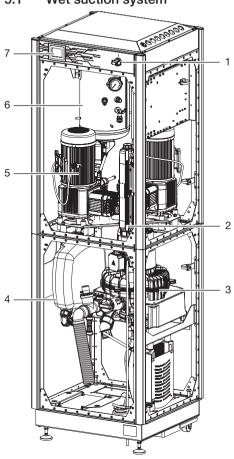

- 1 Main power switch
- 2 Drying unit
- 3 Suction unit
- 4 Pressure equalisation tank
- 5 Compressor unit
- 6 Pressure tank
- 7 Touch display

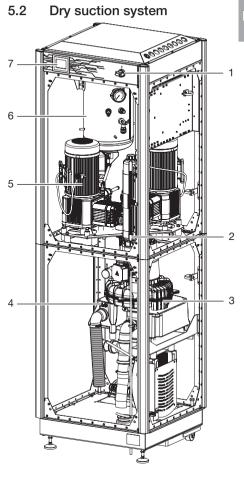

- 1 Main power switch
- 2 Drying unit
- 3 Suction unit
- 4 Condensation separator
- 5 Compressor unit
- 6 Pressure tank
- 7 Touch display



5.3 Functional description

The Power Tower View combines compressed air generation, suction and amalgam separation in one housing. The housing is insulated, ensuring an appropriate level of soundproofing. The unit is operated via a touch display and indicates the various parameters of the individual components.

5.4 Touch display



- Touch display
- 2 Main power switch

The touch display shows unit data, operating statuses and messages, and user-dependent unit settings can be configured.



Assembly



Only qualified specialists or employees trained by Dürr Dental are permitted to install, connect and start using the unit.

6 Requirements



Further information can be found in our suction planning information leaflet.

The room chosen for installation must satisfy the following requirements:

- Closed, dry, well-ventilated room
- It should not be a room made for another purpose (e.g. boiler room or wet cell).
- Set up the unit on a clean, level and sufficiently stable surface (take the weight of the unit into account).
- Set up the units so that the type plate can be easily read and the unit is easily accessible for operation and maintenance.
- Leave sufficient distance from any wall (at least 20 cm).
- Refer to the requirements for environmental conditions in "4 Technical data".

6.1 Pipe materials

Only use vacuum-sealed HT-waste pipes manufactured from the following materials:

- Polypropylene (PP).
- Chlorinated polyvinyl chloride (PVC-C).
- Plasticizer-free polyvinyl chloride (PVC-U),
- Polyethylene (PE).

The following materials must not be used:

- Acrylonitrile-butadiene-styrene (ABS),
- Styrene copolymer blends (e.g. SAN + PVC).

6.2 Hose materials

For waste connections and suction lines only use the following hose types:

- Flexible spiral hoses made of PVC with integrated spiral or equivalent hoses
- Hoses that are resistant to dental disinfectants and chemicals

(j

Plastic hoses will display signs of ageing over time. Therefore, they should be inspected regularly and replaced as necessary.

The following types of hoses must not be

- Rubber hoses
- Hoses made completely of PVC
- Hoses that are not sufficiently flexible

6.3 Installation and routeing of hoses and pipes

- Execute the on-site pipe installation in accordance with the applicable local regulations and standards
- Lay the hose installation of the drains to or from the unit at a sufficient incline.



If incorrectly laid, the hoses can become blocked with sedimentation.

6.4 Information about electrical connections

- Ensure that the electrical connections to the mains power supply are established in accordance with current valid national and local regulations and standards governing the installation of low voltage units in medical facilities.
- The connection to the mains supply must be a fixed connection that cannot be disconnected without the use of tools. Plug-in connections (power outlet/plug) are not permissible.
- Doserve the current consumption of the devices that are to be connected.

Electrical fusing

Version 400 V, 3N~

LS switch 16 A, characteristic C and D in accordance with IEC 60898.

Version 230 V. 1N~

LS switch 25 A, characteristic C and D in accordance with IEC 60898.



6.5 Information about connecting cables

The diameter of the connections depends on the current consumption, length of line and the ambient temperature of the unit. Information concerning the current consumption can be found in the Technical Data supplied with the particular unit to be connected.

The following table lists the minimum diameters of the connections in relation to the current consumption:

Current consumption of unit [A]	Cross-section [mm ²]
> 10 and < 16	1.5
> 16 and < 25	2.5
> 25 and < 32	4
> 32 and < 40	6
> 40 and < 50	10
> 50 and < 63	16

Control cable

24 V protective low voltage for:

- Hose manifold
- Place selection valve
- Spittoon valve

Installation type	Line layout (minimum requirements)
Fixed installation	 Shielded sheathed cable (e.g. (N)YM (St)-J)
Flexible	 PVC data cable with shielded cable sheath- ing, as used for tele- communications and IT processing systems (e.g. type LiYCY)
	or
	 Lightweight PVC control cable with shielded cable sheathing



Connect the shielding of the cables in accordance with the regulations.

7 Installation


7.1 Combining devices safely

Take care when connecting units together or to parts of other systems as there is always an element of risk (e.g., due to leakage currents).

- Only connect units when there can be no question of danger to operator or to patient.
- Only connect units when it is safe to do so and when there is no risk of damage or harm to the surroundings.
- If it is not 100% clear from the unit data sheet that such connections can be safely made or if you are in any doubt, always get a suitably qualified person (e.g. the manufacturer) to verify that the setup is safe.


7.2 Setting up the suction module

- > Set up the lower part of the cabinet.
- Remove the transport protection.
- Attach the cover at the rear.

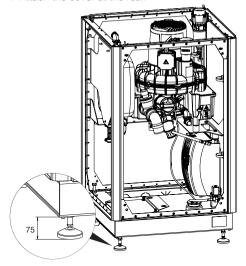


7.3 Install the compressed air

Mount the compressed air module on the suction module using four M5x20 thread-cutting screws. Screw two screws into the left-hand frame and two into the right-hand frame.

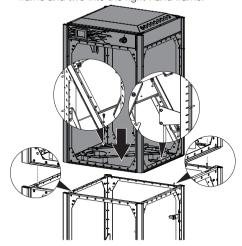

Push on and secure the noise insulation on the suction manifold of the compressor unit.

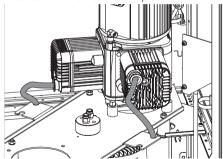

7.5 Positioning the unit

Place the foam over the connections in the floor.

7.4 Installing the compressor unit




CAUTION

Functional disorder in case of incorrect installation

Unit does not start up.


When retrofitting a compressor unit in the case of 400 V variants, the compressor unit must have a star connection.


Position the compressor unit in the unit.

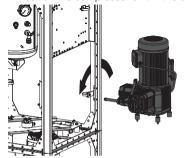


Mount the compressor unit on the retaining plate using four nuts. ΕN

Position the unit over the foam in such a way that the connections are easily accessible.

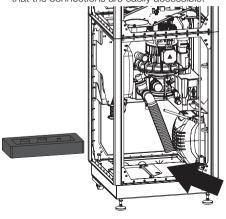


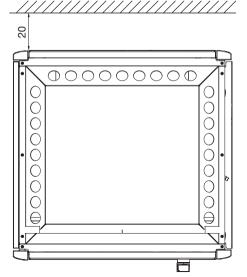

> Align the unit in its final position.

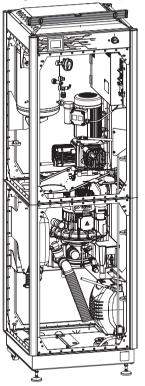

Connect the exhaust air hose and suction hose to the floor connection.

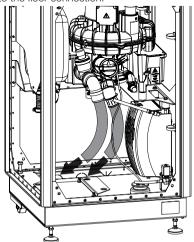

7.6 Installing a second compressor unit (optional)




CAUTION

Functional disorder in case of incorrect installation


Unit does not start up.

- When retrofitting a compressor unit in the case of 400 V variants, the compressor unit must have a star connection.
- Mount the second compressor unit and push on and secure the noise insulation on the suction manifold of the compressor unit.

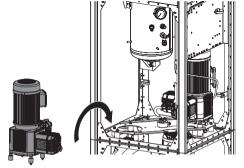


7.7 Install the drying unit

> Install the drying unit.

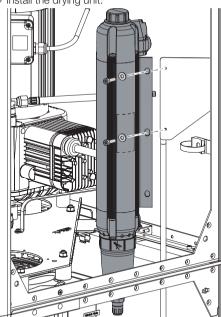


Insert nozzle, see installation instructions 9000-610-49/01. Insert the filter and close the drying unit with the cover.

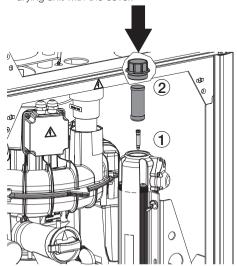




> Connect the drying unit to the compressed air hose coming out of the cooler.

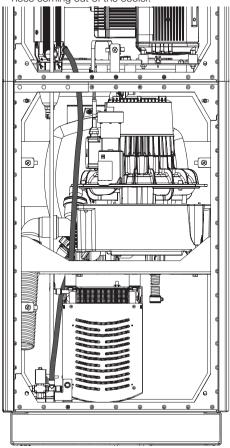


> Connect the drying unit to the pressure tank.

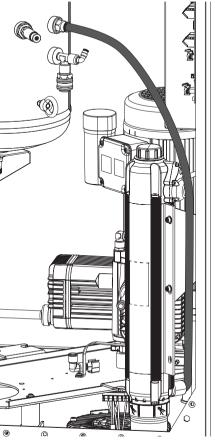


7.8 Additional work for 230 V variant

Operation with two compressor units


> Remove the earth strap.

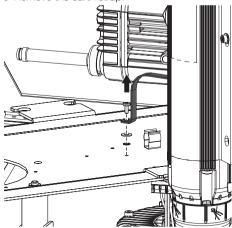


> Remove the clip for L bracket.

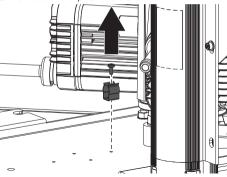


Mount the valve block with earth strap.

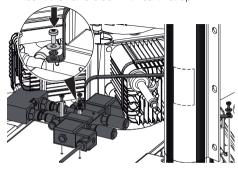

> Establish the compressed air connection between the compressor unit and valve block.

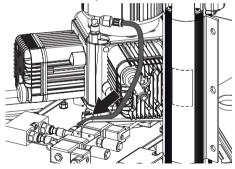


- Mount the second compressor unit and establish the compressed air connection.
- > Route the compressed air hose from the valve block to the cooler and connect it to the cooler.

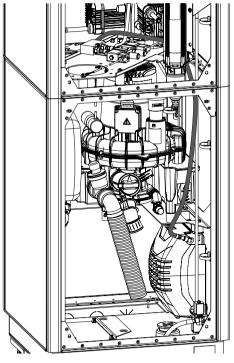


ΕN

> Lay the valve cable from the valve block to the

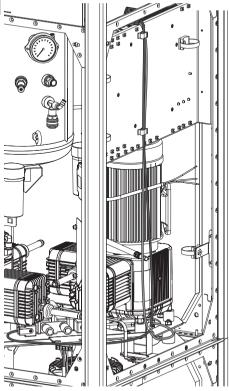


> Establish the electrical connections for the valves (see also "14.2 Version 230 V, 1N~").

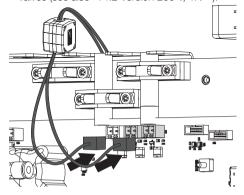


Operation with one compressor unit

Mount the L bracket with the compressed air hose to the cooler.

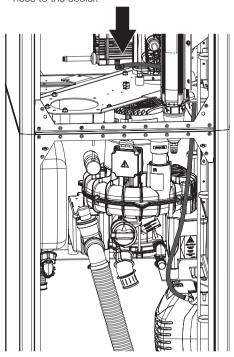


> Establish the compressed air connection for the compressor unit.

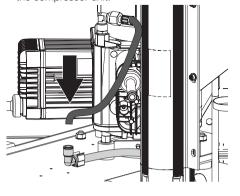


7.9 Install the amalgam separator (optional)


Position the amalgam separator in the lower part of the cabinet.

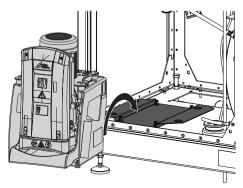


- > Remove the electronic hood.
- Plug two spacers and then BUS adapters into slot X9 of the electronic circuit board.

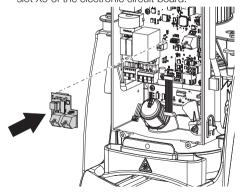

Attach the wastewater connection of the pressure equalisation tank at the inlet of the amalgam separator and the wastewater hose on the amalgam separator outlet.

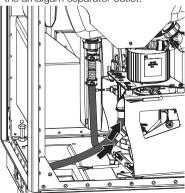

- Position the amalgam separator correctly in the Power Tower View.
- Plug the mains cable into the mains socket on the BUS adapter.

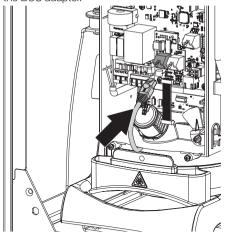


7.10 Connecting the unit

- Make sure that none of the electrical cables leading to the unit are under any mechanical tension.
- Defore start-up, check the mains voltage against the voltage indicated on the type plate (see also "4 Technical data").



- Ensure that the electrical connections to the mains power supply are established in accordance with current valid national and local regulations and standards governing the installation of low voltage units in medical facilities.
- The connection to the mains supply must be a fixed connection that cannot be disconnected without the use of tools. Plug-in connections (power outlet/plug) are not permissible.
- Doserve the current consumption of the devices that are to be connected.
- > Route the prepared cables in the unit, fasten to the pressure tank using Velcro and connect to the control board.

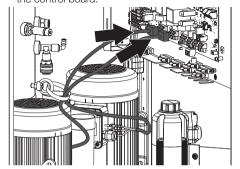

- Mount the cover on the control box.
- Mount the cover in the unit via the floor connections and side panels.

7.11 Two devices in a single compressed air network

In order to connect two units to a single compressed air network,


- the pressure vessels need to be connected to each other;
- the controllers need to be connected to each other.

Connecting pressure vessels

If two devices are connected to a single compressed air network, pressure equalisation must take place between the pressure vessels. To do this, the pressure vessels need to be connected to each other.



So that the pressure can be equalised, no non-return valves must be installed between the pressure vessels.

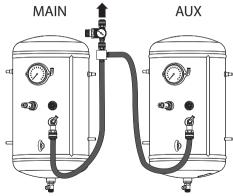

- Main device (MAIN)
- 2 Auxiliary device (AUX)

MAIN / AUX for main device / auxiliary device

The two electronic controllers of the units are connected to each other via a network cable.



When routeing the cables, maintain the correct gaps between control cables and supply cables.

Connect the network cable to the network socket X1 or X5.




8 Commissioning



Perform an electrical safety check prior to commissioning.

- Switch on the main power switch. Compressor will start and switches off at approx. 0.75 MPa.
- In the case of a 230 V variant with one compressor unit, deactivate the second compressor unit in the following menu: System settings > Unit configuration > Compressor > Activate / Deactivate second motor.
- Open the condensate drain valve and drain the condensate.

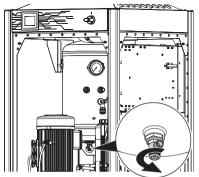

- Switch on the suction machine by removing the suction hose from the hose manifold or by activating the rinsing of the spittoon (spittoon valve, "wet suction system").
- Carry out a functional inspection and check the connections for leaks.
- > For alternating current suction machines, check the direction of rotation.



The VS suction machines continue after running for approx. 30 seconds. In the event of simultaneous start signals for the suction unit and the compressor units, the compressor units start up with a time delay to the suction unit.

Switch off the suction machine by removing the suction hose from the hose manifold or by deactivating the rinsing of the spittoon.



An adequate ventilation by fans in the unit is only possible when all covers are attached. Do not cover the air inlets on the housing.

> Place all covers on the unit.

8.1 Monitoring the unit via the network

The following requirements must be met in order to monitor the unit on the computer:

- Unit connected to the network
- Current monitoring software installed on the computer

Combining devices safely

- The overall safety of the unit and its main performance features are independent of the network. The device is designed for operation independent of a network. However, some of the functions are not available in this case.
- Incorrect manual configuration can lead to significant network problems. The expert knowledge of a network administrator is required for configuration.
- The data connection utilises part of the bandwidth of the network. Interactions with other medical devices cannot be completely ruled out. Apply the IEC 80001-1 standard for risk assessment.
- The device is not suitable for direct connection to the public Internet.

Network configuration

Various options are available for network configuration:

- ✓ Automatic configuration via DHCP (recommended).
- ✓ Automatic configuration via Auto-IP for direct connection of unit and computer.
- ✓ Manual configuration.
- > Configure the network settings of the unit using the software or, if available, the touch screen.
- Check the firewall and release the ports, if applicable.

Network protocols and ports

Port	Purpose	Service
45123 UDP, 45124 UDP	Unit recognition and configuration	
1900 UDP	Service indicator	SSDP/ UPnP
502 TCP	Device data	
514 ¹⁾ UDP	Event log data	Syslog
22 TCP, 23 TCP	Diagnosis	Telnet, SSH
123 UDP	Time	NTP



Port Purpose Service

1) The port may vary depending on the configuration.



9 Operating the touch screen



NOTICE

Damage to the touch screen due to incorrect handling

- Only touch the touch screen with your fingertips.
- Do not use a sharp instrument (e.g. ballpoint pen) to operate the touch screen.
- Protect the touch screen against water.
- Operate the touch screen by tapping it with a fingertip to select a button or input field.



9.1 Navigating

If the contents of the window cannot be completely displayed on the touch screen, a scroll bar appears.



Tap or to move the displayed section of the window.

9.2 Using menus


Use the buttons to switch to other menus.

- Tap \equiv to switch to **Settings**.
- > Tap Next to scroll.
- ➤ Tap ♠ to switch to the start screen.
- ➤ Tap ← to switch to the previous menu level.

Tap ② to switch to the detail view for the alarm texts.

9.3 Calling up messages on the touch screen

Messages are divided into the following categories:

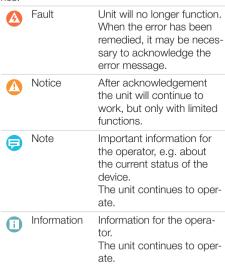




10 Operation

10.1 Switch on the unit.

> Switch on the unit at the main power switch.

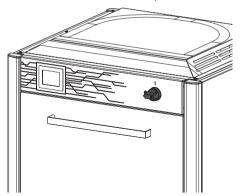

10.2 Device information

Use the menus *Unit information > Unit data* and *Unit information > Unit usage data* to call up information including:

- IP address
- Serial number
- Firmware version

10.3 Configuring the units

Use the *Parameters* menu to configure settings including:

- Compressor pressure range
- Suction unit lag time

Use the *System settings* menu to configure settings including:

- IP address
- Language

10.4 Amalgam separator

> Tap ŵ

The fill level of the amalgam collecting container and the pressure in the pressure tank are displayed.

> Tap Next

The humidity and the temperature of the electronics are displayed.

ΕN



11 Maintenance



Prior to working on the unit or in case of danger, disconnect it from the mains.

11.1 Compressor unit and drying unit



WARNING

Risk of infection due to burst filters

Particles enter the compressed air network and can therefore enter the mouth of the patient.

> Replace filters in accordance with the maintenance schedule.



Wear protective equipment to avoid any risk of infection (e.g. liquid-tight protective gloves, protective goggles, face mask).

Maintenance schedule



NOTICE

Risk of damage to the unit due to blocked filters


Continuous running due to reduced delivery. Damage to the unit due to burst filters.

> Replace filters in accordance with the maintenance schedule.

Maintenance interval	Maintenance work
Annually	 Replace the air intake filter in the compressor unit – do this every six months if there is a high concentration of dust. Replace the fine or virus bacteria filter. Replace the sintered filter.
In accordance with national law	 Check the safety valve. Carry out recurring safety inspections (e.g. pressure tank inspections, electrical safety inspections) in accordance with applicable national laws.

Replacing the air intake filter

- Switch off the device.
- Pull off the noise reducer from the air intake filter
- Remove the air intake filter.
- Insert a new air intake filter.
- Push on the noise reducer onto the air intake filter.

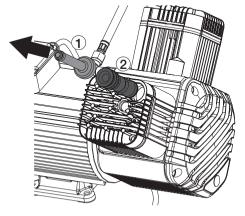

Replacing the filter of the drying unit

Fine or virus bacteria filter

- > Unscrew and remove the filter cover.
- Remove the filter.
- Insert a new filter.
- > Replace the filter cover and close.



Sintered filter

- > Unscrew and remove the filter housing.
- > Remove the filter.

- Insert a new filter.
- > Replace the filter housing and close.



Wear parts and replacement parts

The following working parts must be replaced at regular intervals (refer also to "Maintenance"); these articles do not bear the CE mark;

Air intake filter	0832-982-00
Fine filter	1610-121-00
Virus bacteria filter	1650100172
Sintered filter	1650-101-00
Coalescence filter	1650200323



To configure the required filters or filter sets, you can also use our filter configurator at:

www.duerrdental.com/filterkonfigurator




Any repairs exceeding routine maintenance may only be carried out by qualified personnel or our service.



Information about replacement parts is available from the portal for authorised specialist dealers at:

www.duerrdental.net

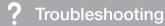

ΕN

12 Tips for operators and service technicians



Any repairs exceeding routine maintenance may only be carried out by qualified personnel or our service.




Prior to working on the unit or in case of danger, disconnect it from the mains.

Error	Possible cause	Remedy
Direction of rotating field incorrect	The phases of the supply voltage cable have been mixed up, the suction unit motor may be rotating in the wrong direction	Switch the phases of the sup- ply voltage cable.


Appendix

13 Connection media plan

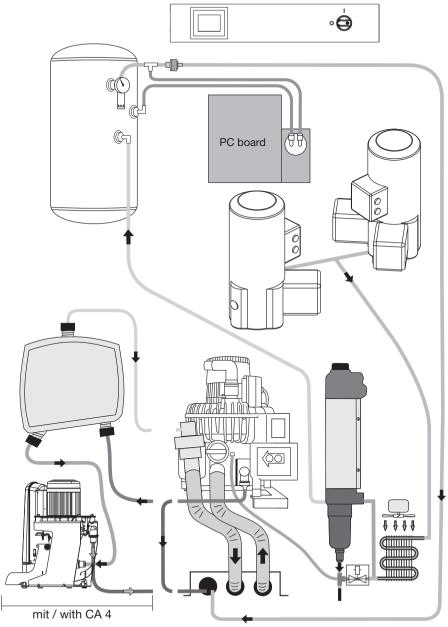


14 Electrical connection plan

14.1 Version 400 V. 3N~

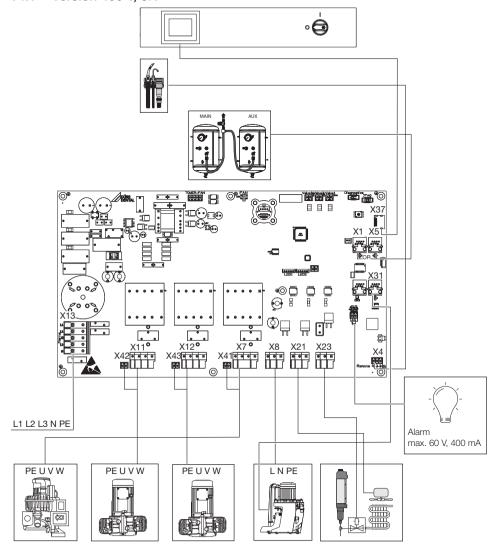

ΕN

Version 230 V, 1N~

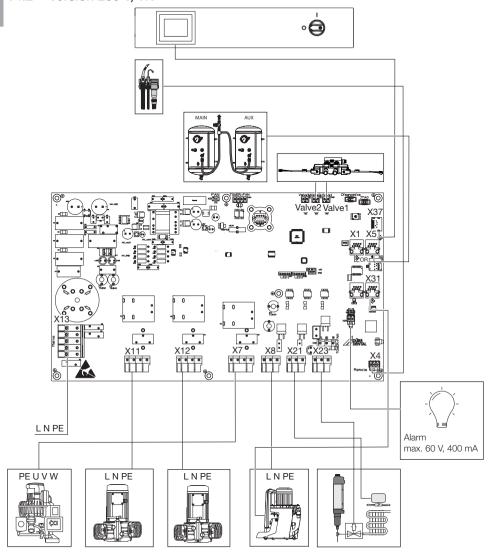

15 Menu structure

15.1 Settings

10.1 Octungs			
Level 1	Level 2	Level 3	Level 4
Access level ¹	Operator		
	Administrator		
	Service Technician (PIN)		
Device Information ¹	Device data		
	Device usage data		
System settings ¹	Language 1	German (DE)	
		English (EN)	
	Date / time ²	Time zone	
		Date	
		Time	
	Network ²	DHCP	
		IP address	
		Netmask	
		Gateway	
	System of units ²	Metric	
		Imperial	
	Unit configuration ³	Amalgam separator ³	Remove amalgam separator
		Compressor ³	Activate / Deactivate second motor
		Cluster settings ³	Auxiliary unit
			Main unit
	Factory settings ³	Clear message history	
Parameters ³	Start-up pressure ³		
	Cut off pressure ³		
	Lag time		
Message history ²			
Maintenance ³	Maintenance completed		
Service menu ²	Compressor	Motors	Start motors
		Emergency mode ²	Reset
	Amalgam separator ³	Start motor	
	Suction unit ³	Start motor	
	Oddion drift		





- Visible from access level Operator or higher
- Visible from access level *Administrator* or higher
- ³ Visible from access level *Service Technician* or higher



The access level "Service Technician" offers additional functions that must only be adjusted by Dürr Dental or by individuals/companies authorised to do so by Dürr Dental. Corresponding documents for this can be downloaded from www.duerrdental.net.

This document confirms that a qualified handover of the medical device has taken place and that appropriate instructions have been provided for it. This must be carried out by a qualified adviser for the medical device, who will instruct you in the proper handling and operation of the medical device.

Product name	Order number (REF)	Serial number (SN)		
 □ Visual inspection of the packaging for any damage □ Unpacking the medical device and checking for damage □ Confirmation of the completeness of the delivery □ Instruction in the proper handling and operation of the medical device based on the operating instructions 				
Notes:				
Name of person receiving instru	uction: Signature:			
Name and address of the qualified adviser for the medical device:				



Date of handover:	Signature of the qualified adviser for the medical device:



Hersteller / Manufacturer:

DÜRR DENTAL SE Höpfigheimer Str. 17 74321 Bietigheim-Bissingen Germany Fon: +49 7142 705-0

Fon: +49 7142 705-0 www.duerrdental.com info@duerrdental.com

